


# LED Pet Piezo™ Scaler

The **LED PET PIEZO™ SCALER** is a compact piezoelectric scaler. The unit features Auto Gain Control which stabilizes and adjusts power settings for optimal performance. The **LED PET PIEZO™ SCALER** also has a large power dial surrounded by LED lit indicators for readability and easy access. The handpiece has a bright LED light that illuminates the work area for easy use while scaling.

**FEATURES OF THE LED PET PIEZO™ SCALER:**

- LED Scaler Handpiece
- Designed to accommodate KLAW™ Tips
- Includes 4 scaler tips
- External Water Source: Inline or Pump

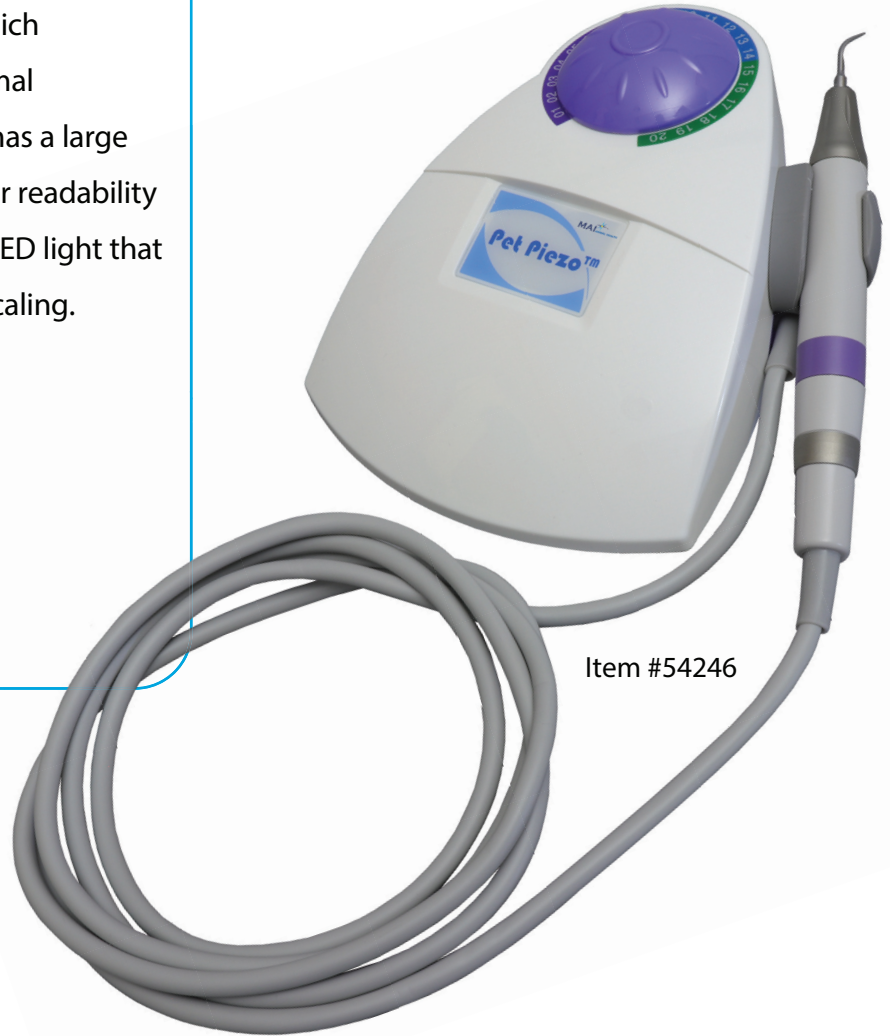

Item #54246

Item	Product Description	UOM
54246	LED Pet Piezo™ Scaler	Each

**STRONG PERFORMANCE,  
DEPENDABLE QUALITY**

KLAW™ Tips are ultrasonic instruments ideal for cleaning deep developmental grooves in cat canine teeth and dog premolar teeth. These instruments are specifically designed to cause less tooth damage when cleaning and to increase patient comfort during dental cleanings. The range of KLAW™ Tips were created exclusively for use with the LED Pet Piezo™ Scaler. The tips are coated with Zirconium Nitrate for durability and have ported tips to ensure water is being applied to the area being cleaned. KLAW™ Tips are small and easily changeable onto the piezoelectric scaler handpiece.

**KLAW™ 1:**

A universal small instrument with a curved scaler tip designed for root surfaces. This curette design will eliminate hand scaling repetitive motion injuries and is useful for subgingival soft tissue curettage.

Item	Product Description	UOM
54830	KLAW™ 1 - Piezo Tip	Each



Item #54830

**KLAW™ 2:**

A small instrument with a left 30° angle near the tip for ease of scaling difficult to reach areas. Designed for right mandibular buccal surfaces, left mandibular lingual surfaces, left maxillary buccal surfaces and left maxillary lingual surfaces.

Item	Product Description	UOM
54831	KLAW™ 2 - Piezo Tip	Each



Item #54831

**KLAW™ 3:**

A version of KLAW™ 2 with a right 30° angle near the tip for ease of scaling difficult to reach areas. Designed for left mandibular buccal surfaces, left mandibular lingual surfaces, left maxillary buccal surfaces and left maxillary lingual surfaces.

Item	Product Description	UOM
54832	KLAW™ 3 - Piezo Tip	Each



Item #54832

**KLAW™ 4:**

A straight universal instrument for use on all tooth surfaces which is tapered with a tip angle of 110°. Ideal for interdental areas with a sharpened blade useful for root scaling.

Item	Product Description	UOM
54833	KLAW™ 4 - Piezo Tip	Each



Item #54833

**KLAW™ 6:**

A straight universal instrument for use on all tooth surfaces which is tapered with a tip angle of 110°. Ideal for interdental areas with a sharpened blade useful for root scaling.

Item	Product Description	UOM
54334	KLAW™ 6 - Piezo Tip	Each



Item #54834