

# Olympic 3-in-1 Treatment Table™

A complete solution for all your weighing,  
treatment and transport needs.

**NEW**

## All-Purpose Table

- Transport
- Treatment
- Exam



Scan to watch video  
demonstration

or visit us at [www.olyvet.com](http://www.olyvet.com)



### Built-in Precision Scale – always at the ready

The Olympic 3-in-1 Treatment Station includes a precision scale that swivels and tilts so it is always at an optimal viewing angle. The self-calibrating scale maintains a  $\pm 0.1$  lb accuracy – whether you are weighing a 250 lb pot belly pig or a 1 lb bird. And the weight averaging function locks in quickly, even on active patients. Also includes a 'last weight' recall button.

### Customizes in seconds with 6 drop-in sockets

This work station is designed for efficiency and ergonomic comfort, standing or seated. All your instruments, monitors and other equipment are at your fingertips to save you time. And since all the accessories are mounted on the table, it eliminates floor clutter and saves space. Create the optimal workspace for each procedure with an Olympic IV Pole, Laptop Tray, Swiveling Instrument Tray, LED Procedure Light or Equipment Shelf.

### Reduces handling of compromised patients

The Olympic Treatment Station is designed to navigate easily from exam to treatment, x-ray, prep, surgery or recovery – and even outdoors. By combining patient treatment and transport, it reduces the handling of injured or compromised patients and saves time. The built-in scale is especially helpful when medications need to be administered and there is no time to use the floor scale.

### Patient-friendly, non-slip Corian® tabletop

Table features an optional Corian® top, the ideal surface for modern hospitals. Corian feels warmer than stainless, is non-skid and non-reflective so patients are more relaxed during exams and treatment. Features anti-drip edge around perimeter.

### Convenience Features

The 3-in-1 Treatment Station features an on-board power strip, 8 tie-down cleats, push handles, 2 foot switches and a long-life rechargeable battery.

## Precision scale is built-in



### Equipment Shelf

For monitors and other equipment. Easy to see and reach. Mounts at head of station. Electrical equipment connects directly to power outlets on table (see below).



### Electrical Outlets

7 electric outlets for connecting all your accessories and equipment. Conveniently located at end of table near equipment shelf.



### Laptop

Your laptop is securely held on a special swing-arm so it's just where you want it. Can be quickly moved to any drop-in socket.



### Instrument Tray

On swing-arm. Fits sockets at 4 corners. Instruments always at your fingertips, raises and lowers with the table. 13" x 10".



### IV Pole

Drops into sockets at corners of station. Eliminates floor stand. Travels with patient. Can support fluid warmer.




### Procedure Light

Daylight color LED. Adjustable brightness and focus. Connects to electrical outlets on the 3-in-1 Treatment Table.




### Scale

Scale under Corian® top. Large, easy to read. Swivel mounted. Can weigh active animals. Accuracy +/- 0.1 lb. Lifetime quality.

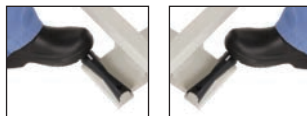

### 2 Foot Switches

For faster and convenient height adjustment, there are foot switches on both sides of the table.

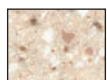

**Sahara**  
51680



**Sandstone**  
51681



**Platinum**  
51682

### Corian® non-slip surface



Table top surface is lifetime Corian® from DuPont, the ideal surface for modern hospitals. Feels warmer than stainless steel, is non-skid so your patients remain calm, wipes clean with disinfectants and features a drip edge. Choose from 3 standard colors/patterns.

## For all procedures and transport to or from the surgical suite



### Ergonomic and efficient

The 3-in-1 Treatment Table offers ergonomic comfort with full-range height adjustment. The maximum height is an extra tall 40". Work standing or seated. Improved patient presentation. The rechargeable electric lift is quiet and smooth. Two foot switches, one on each side.

### Patient transport – multi-purpose

An easy-roll patient transport. It has steering handles and is light weight so it turns easily in tight spaces. Now you can seamlessly move patients between treatment, x-ray, prep, surgery, recovery.



### Specifications

**Dimensions** – Overall 23" x 58",  
Tabletop 23" x 53"

**Height range** – 13" to 40"

**Lift capacity** – 250 lbs

**Casters** – Large 5", full swivel, locking

**Tie-down cleats** – Under both  
sides of table

**Net weight** without accessories,  
95 lbs

**Electrical** – 120V/60Hz

## Order From Your Distributor

### Models with Corian® Table Top

**Olympic 3-in-1 Treatment Table - Complete** ..... Model 51601  
Includes Scale, Corian® Top, IV Pole, Swiveling Instrument Tray, Equipment Shelf, Ultra-Bright LED Procedure Light and Swiveling Laptop/Monitor Tray.

**Olympic 3-in-1 Treatment Table - Base Model with Scale...** Model 51602  
Includes Scale, Corian® Top and 4 sockets ready for accessories.

**Olympic 3-in-1 Treatment Table - Base Model without Scale...** Model 51603  
Includes Corian® Top and 4 sockets ready for accessories.

### Models with Stainless Steel Table Top

**Olympic 3-in-1 Treatment Table - with Scale** ..... Model 51604  
Includes Scale, Stainless Steel Top and 4 sockets ready for accessories.

**Olympic 3-in-1 Treatment Table - without Scale** ..... Model 51605  
Includes Stainless Steel Top and 4 sockets ready for accessories.

### Accessories

**IV Pole** ..... Model 51611

**Swiveling Instrument Tray, 13" x 10"** ..... Model 51612

**Equipment Shelf, 11" x 20"** ..... Model 51613

**Ultra-Bright LED Procedure Light** ..... Model 50970

**Advanced LED Procedure Light** ..... Model 50971

**Swiveling Laptop/Monitor Tray** ..... Model 51615

**3-year Olympic Warranty. All parts and labor.**