

Olympic Hi-Lo Wet Table™

Better presentation – More comfortable
The Only Wet Table With Adjustable Height And Swing Arms



Wet Table Features

This is the modern wet table that replaces wooden built-ins. Offers height adjustment and useful accessories, also converts easily to a fully equipped dental table.

Enclosed Electric Lift – Precision height adjustment. Capacity 250 lbs. Clean, fully enclosed.

Ergonomic Comfort – Work seated or standing on large or small patients. Meets ergonomic standards.

Better Patient Presentation – With adjustable height. More efficient, faster treatments.

Large Surface – 60"x 24", 2 removable grates.

Rock Solid – Reinforced all welded steel construction. No bounce or wobble. A better work platform.

Leveling Feet – Assure solid work surface.

Swing Arms – Choice of one or two Swing Arms with drop-in sockets for IV Pole, Light, etc. Hold 60 lbs each. Go up/down with table. Arms lift in/out of brackets on column without tools.

Choice Of Accessories – The only wet table that provides useful accessories for both wet and dry treatments and dental procedures.


Install Anywhere – Fits any standard plumbing connection. Unique design – no noisy pump, no external hoses. Installs near wall or free standing. Not a built-in, can be moved in future.

Saves Space and Improves Efficiency



Scan to watch video demonstration

or visit us at www.olyvet.com

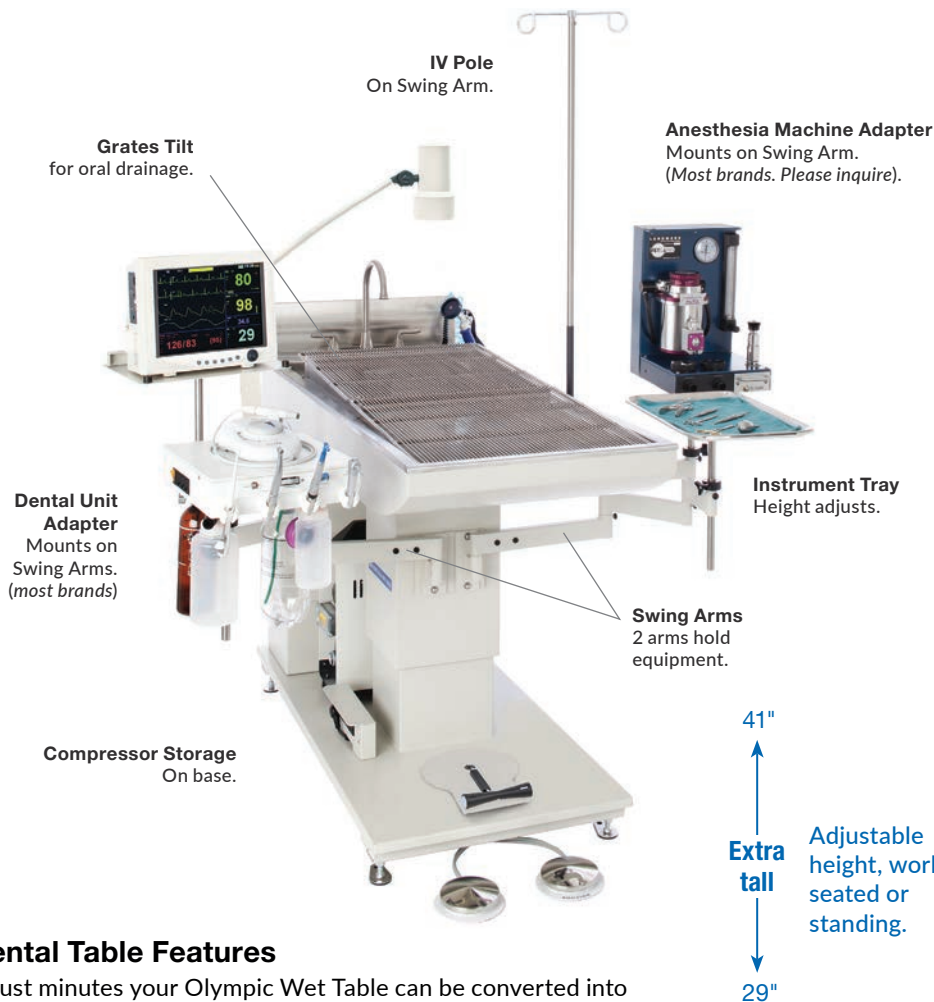

LED Procedure Light
Bright, focusing,
flexible.



Laptop Platform
Your laptop
held securely.



4 Electrical Outlets



Grate Tilt
Fits under grates.

Dental Table Features

In just minutes your Olympic Wet Table can be converted into a fully equipped dental table. All your instruments and equipment are positioned at your fingertips, ready for faster, more efficient procedures.

Ergonomic Height Adjustment – You can work more comfortably seated or standing. Less tiring for long dental procedures.

Improved Presentation – Another benefit of height adjustment. Better visualization and access to the oral site.

Use Your Current And Future Equipment – This table is designed as a universal platform. You can use almost any dental unit or anesthesia machine.

Accessories – A choice of several helpful dental accessories currently available (next column).

Solid Work Surface – The heavy-duty construction of the Olympic Hi-Lo Wet Table provides a rigid platform for dental procedures.

Reduces Floor Clutter – Mounting equipment on the table saves floor space and congestion.

Swing Arms – For dental use, two Swing Arms are recommended. These articulating arms each have 2 drop-in sockets for your dental equipment and accessories. Each arm can support 60 lbs. The arms lift in/out of sockets on the lift column and can be quickly changed without tools.

2 Tables In One – The Olympic Hi-Lo Wet Table not only provides a superior wet-dry treatment table but also a fully equipped and efficient dental table.

Order From Your Distributor

Olympic Hi-Lo Wet Table - Complete System..... Model 51003

Complete for treatments or dentistry. Includes 2 Swing Arms with 2 sockets each, height-adjustable Instrument Tray, IV Pole, Ultra-Bright LED Procedure Light, Laptop/Monitor Tray, Grates, Grate Tilt Insert and Stainless Steel Hard Top, all plumbing fixtures. Overall 24" x 60" x 41". Weight 360 lbs. 120V/60Hz.

Olympic Hi-Lo Wet Table - Base Unit Model 51001

Includes 2 Grates and all plumbing fixtures. Includes 2 sockets for mounting Swing Arms and accessories.

Accessories

Accessories quickly changed for wet table or dental use

IV Pole Model 51081

Swing Arm with 2 sockets Model 51086

Special Use Swing Arm Model 50986

For mounting some brands of large anesthesia units

Swiveling Laptop/Monitor Tray Model 50988

Instrument Tray, 13" x 10", height-adjustable Model 51087

Instrument Tray, 19" x 13", height-adjustable Model 51884

Grate Tilt Model 51082

Stainless Steel Dry Top Model 51085

12" Arm Extension (adds 3rd socket) Model 50996

Ultra-Bright LED Procedure Light Model 50970

Advanced LED Procedure Light Model 50971

Adapters available to mount most anesthesia and dental units.

3-year Olympic Warranty. All parts and labor.